



A COMMON INTRODUCTION

McStas & McXtrace

Monte Carlo ray-tracing for neutron and X-ray instrumentation

McStas — neutrons

McXtrace — X-rays



PaNRAID School · supported by PEPR DIADEM



Who's behind this school



PaNRAID is organised with support from the France 2030 PEPR DIADEM programme, training researchers to combine McStas and McXtrace ray-tracing simulations with AI methods.

PARTNERS & SUPPORTING ORGANISATIONS



What we'll cover

1

Shared origins & philosophy

One codebase family, two physics engines

3

Anatomy of a simulation

Instruments, components, and the ray state

5

Neutrons vs. X-rays

Where the two probes really differ

2

The Monte Carlo ray-tracing idea

Why rays, weights, and randomness

4

Under the hood

DSL, code generation, and coordinate placement

6

Community, use & getting started

License, ecosystem, and next steps



One family, two probes

McXtrace was built on proven McStas technology — same architecture, adapted for X-rays.

Commonality: same domain-specific language, same component-based architecture, same open, community-driven development model.

1998

2000s

2009

Today

McStas is born

Grows into a general tool

McXtrace forks off


Parallel development

K. Nielsen & K. Lefmann, Risø.
Originally for triple-axis
spectrometer simulation.

Adopted by DTU Physics, ILL,
PSI, University of Copenhagen,
ESS DMSC. GPL open source.

“Powered by McStas
technology.” Re-targeted at
X-ray optics & beamlines.

Shared kernel philosophy,
separate physics & component
libraries, active communities.

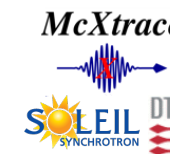




Monte Carlo ray-tracing, in a nutshell

A statistical method born at Los Alamos in the 1940s (*von Neumann, Ulam, Metropolis*) for solving particle-transport problems by numerical experiment.


- Both codes trace many independent rays (neutrons or photons) one at a time through the instrument.
- Propagation in free space is deterministic — straight lines (plus gravity, for slow neutrons).
- Where physics is probabilistic — scattering, absorption, reflection — outcomes are Monte Carlo sampled.
- Each ray carries a statistical weight rather than representing one real particle, so rare but important events can be sampled efficiently.

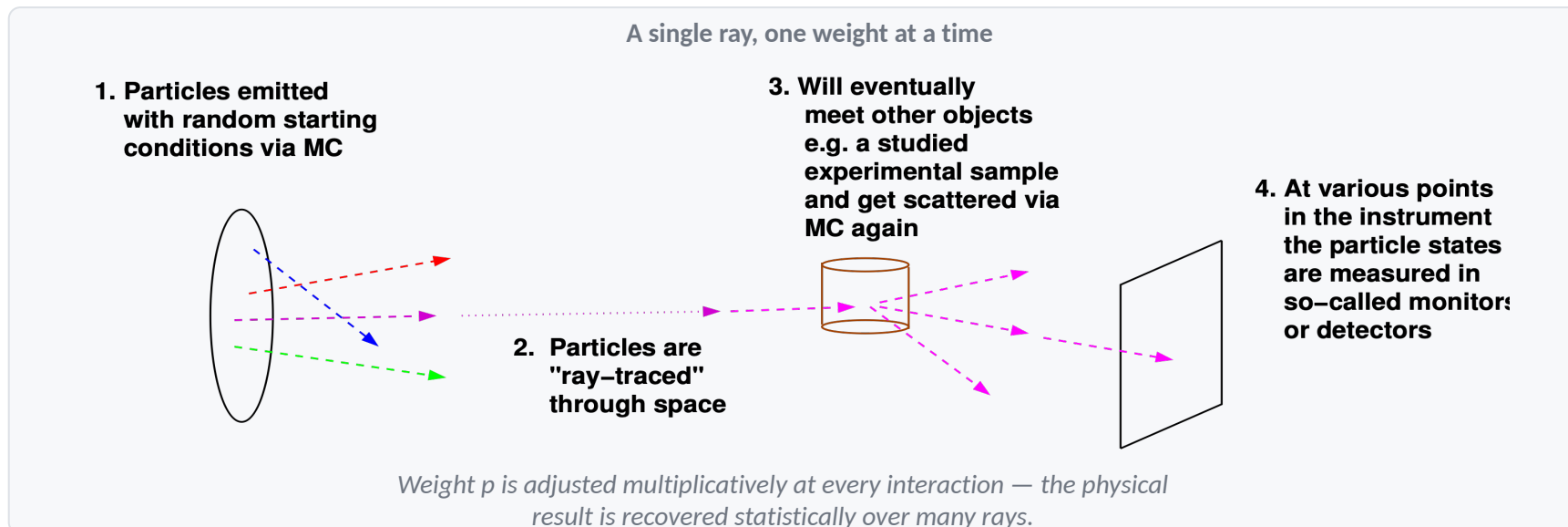



Anatomy of a simulation

An instrument is an ordered chain of components in a shared lab coordinate system. Rays flow through it, changing state at each stop.

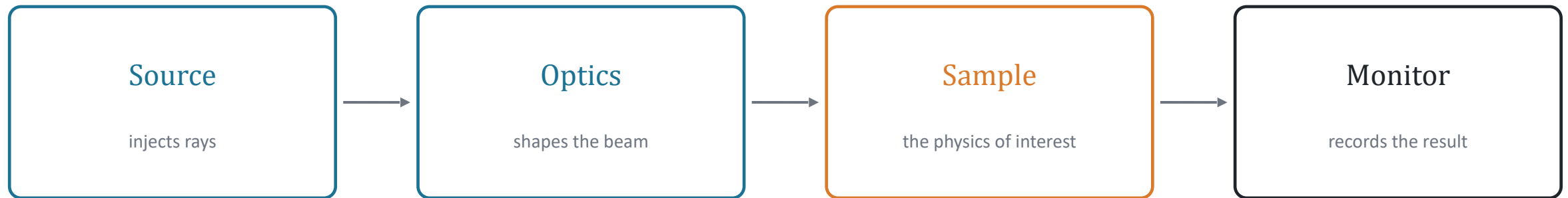

Instrument: defines the lab coordinate system. **Components:** devices/features available in the instrument, each with its own function. **Rays:** passed from one component to the next, changing state along the way.

This chain — and the DSL that describes it — is identical in McStas and McXtrace.

What travels through the beamline: the ray

Neutron ray · McStas



r position (x, y, z)


v velocity vector

t time

S spin / polarisation

p statistical weight

Photon ray · McXtrace

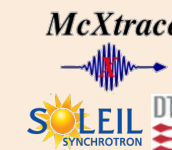

r position (x, y, z)

$\mathbf{k}$ wave vector

ϕ, E phase & polarisation (E-field)

t time

p statistical weight

Same idea, different carrier: both rays are position + time + statistical weight. The neutron adds mass, velocity and spin; the photon adds a wave vector, phase and field polarisation.



Component classes

Every instrument, in either code, is built from the same four component classes.

Class	McStas (neutrons)	McXtrace (X-rays)
Sources	Reactor / spallation moderator	Undulator, wiggler, bending magnet
Optics	Guides, mirrors, choppers, collimators	Mirrors, monochromator crystals, multilayers, slits
Samples	Powders, crystals, liquids, magnetism	Powders, crystals, multilayers, absorbing phantoms
Monitors	PSD / histograms / event lists	PSD / histograms / event lists (same concept)


Common to all components: they set, manipulate, or measure the state of the ray — and each is written as a self-contained, reusable code module.




From instrument file to running simulation



Same toolchain in both codes — the DSL grammar and code generator are shared.

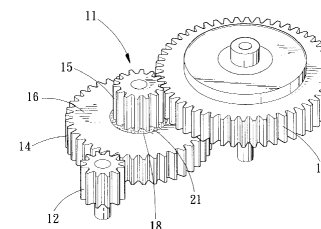

Three levels of code, three audiences

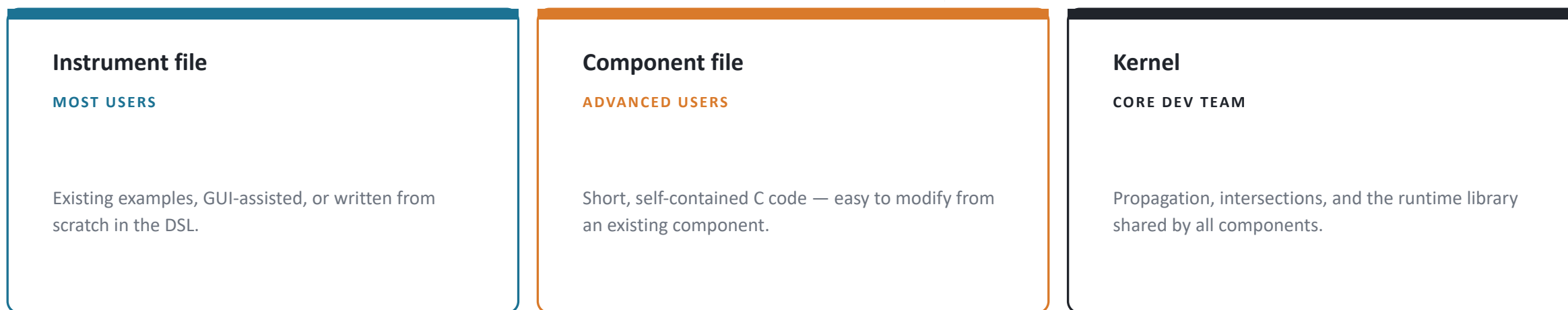





Placing components in space

Every component is placed with the same two statements, in both codes:

```
AT (x, y, z) RELATIVE <ref>  
ROTATED (rx, ry, rz) RELATIVE <ref>
```

- <ref> can be ABSOLUTE (the lab frame), another component's name, or PREVIOUS.
- z runs along the beam, y is vertical, x completes a right-handed frame.
- Component authors choose whatever local origin is physically natural — the kernel handles the coordinate transforms between them.
- AT and ROTATED can be fully decoupled, referring to different RELATIVE <ref>

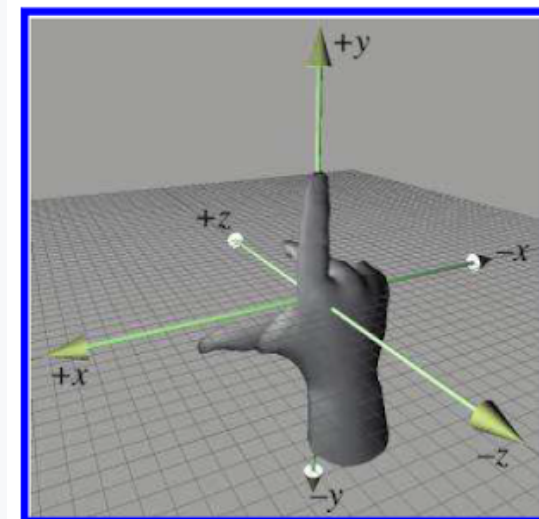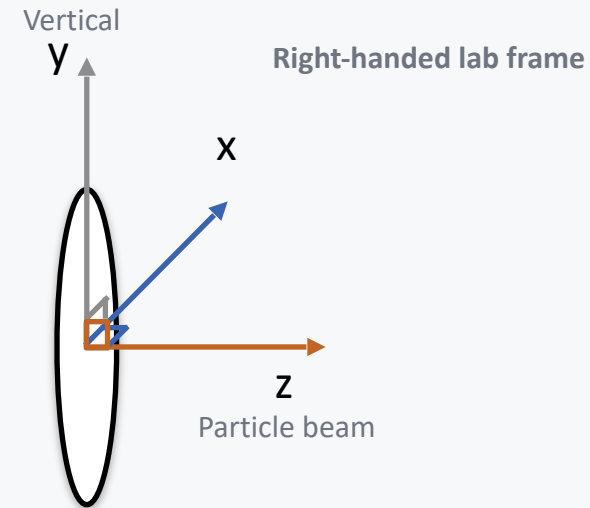

Neutrons vs. X-rays: key differences at a glance

	Neutrons · McStas	X-rays · McXtrace
Probe	Neutral particle — mass, spin	Massless electromagnetic wave
Typical facility	Reactor or spallation source	Synchrotron or free-electron laser
Energy scale	meV (thermal / cold)	keV (hard X-ray)
Dominant interaction	Nuclear & magnetic scattering; deep penetration	Electronic scattering; surface-sensitive, strongly absorbed

What stays the same: the ray-tracing method, the DSL, the component architecture, and the statistical-weight formalism are identical — only the physics inside the components differs.



Open, portable, community-built

GPL

Open-source license
for both packages

3

Platforms: Linux,
macOS, Windows

175+

Components in the
shared-style libraries


30–40%

Of new components
come from users

- Runs everywhere from a laptop to HPC clusters and GPUs (OpenACC) — the same generated C code, recompiled per target.
- New components automatically fit the framework once written in structured, standard C.
- Contributions require at least one test instrument and a documented header — then they join the shared contrib/ library.
- Permanent staff at DTU Physics maintain the core; a wider network of institutes and students keeps both libraries growing.

Who builds & maintains them

A shared core team at DTU Physics, plus partner institutes around each package.

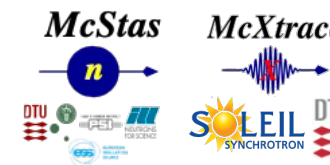


McStas partner network



- DTU Physics (core development)
- European Spallation Source (ESS DMSC)
- University of Copenhagen
- Institut Laue-Langevin (ILL), Grenoble
- Paul Scherrer Institute (PSI)

McXtrace partner network

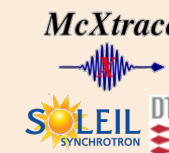

- Synchrotron SOLEIL
- DTU Physics (core development)
- DTU Space
- University of Copenhagen

Further both codes benefit from the greater network of users across a plethora of use-cases and applications.

Several developers work across both projects — shared infrastructure means fixes and features often benefit both codes at once.



What are they used for?

1

Instrument design

Planning and optimising new beamlines before they're built

2

Construction support

Verifying components and layouts during commissioning

3

Virtual experiments

Letting prospective users try an experiment before travelling to a facility

4

Data analysis

Comparing measured data against simulated instrument response

5

Teaching

Training students and new users in instrumentation & scattering methods

Both communities use their code across the full lifecycle — from first sketch to published data.



A shared, interoperable ecosystem

1

MCPL

The Monte Carlo Particle List format lets McStas, McXtrace and other transport codes (MCNP, Geant4, PHITS, ...) exchange rays and couple simulations together.

2

Graphical interface

A common Python-based GUI supports building, running and browsing results for instruments in either code.

3

Plotting tools

A shared plotting toolchain (HTML/D3, pyqtgraph, matplotlib backends) visualises monitor output from both packages.

4

Packaging

Conda-forge distributions and cross-platform installers keep both codes easy to install on Linux, macOS and Windows.

Shared interfaces, tools and utilities



The image displays several screenshots of the McStas and McXtrace software interfaces. The top-left window, titled 'mcgui-py', shows a 3D visualization of a neutron beam passing through a series of mirrors and lenses. The bottom-left window, titled 'templateSANS2_20230518_105746', shows a 3D visualization of a neutron beam passing through a series of mirrors and lenses. The right side of the image shows a 'templateSANS2_generated' window with a Python script for configuring the instrument. The bottom-right window shows a 'Log Console' with a log of the simulation results.

McStas code generator version 3.3 (Mar. 30, 2023)
 Copyright (C) DTU Physics and Risø National Laboratory, 1997-2023
 Additions (C) Institut Laue Langevin, 2003-2019
 All rights reserved
 Components are (C) their authors, see component headers.
 - Component library: McStas-3.3, built Mar. 30, 2023

templateSANS2_20230518_105746
 How To: Inst... Forums... scraping.js... eBay: Neue Un... igen Preisen... Lastbiller Ins... ures... Page 9... Wengo Modelb... ing Nielsen... Jernbanen.dk... _ods lastbiller... Danske jernba... _ Sporskiftet... https://prod..._lang_code=DK

templateSANS2_generated
 Python 3 (pykernel)

```
[ ]: # Comp instance source, placement and parameters
source = instr.add_component('source', 'Source_simple', AT=[0, 0, 0], AT_RELATIVE='arm', ROTATED=[0, 0, 0],

[ ]: source.radius = '0.02'
source.yheight = '0'
source.xwidth = '0'
source.dist = '6'
source.focus_xw = '0.01'
source.focus_yh = '0.01'
source.E0 = '0'
source.dE = '0'
source.lambda0 = 'lambda'
source.dlambda = 'dlambda'
source.flux = '1e8'
source.gauss = '0'
source.target_index = '+ 1'

[ ]: # Comp instance coll1, placement and parameters
coll1 = instr.add_component('coll1', 'Slit', AT=[0, 0, 3], AT_RELATIVE='arm', ROTATED=[0, 0, 0],

[ ]: coll1.xmin = 'UNSET'
coll1.xmax = 'UNSET'
coll1.ymin = 'UNSET'
coll1.ymax = 'UNSET'
coll1.radius = '0.005'
coll1.xwidth = 'UNSET'
coll1.yheight = 'UNSET'

[36]: # Overview plot:
ms.make_sub_plot(data)
```

Log Console
 + Add Checkpoint Clear Log Log Level: Simple 0 1 2

Where to go next



<https://github.com/mccode-dev/PaNRAID/tree/main/Software> has a **.yaml** file
with McStas, McXtrace and the PaNRAID selection of AI libs / tools



McStas

WEBSITE

mcstas.org

MAILING LIST

mcstas-users@mcstas.org

INSTALL

conda-forge · cross-platform installers

SOURCE

GitHub — McStasMcXtrace/McCode

McXtrace

WEBSITE

mcxtrace.org

MAILING LIST

mcxtrace-users@mcxtrace.org

INSTALL

conda-forge · cross-platform installers

SOURCE

GitHub — McStasMcXtrace/McCode

Alternative route just for the MC

`micromamba env create -n panraidMC mcstas mcxtrace -y`
(`apt install mcstas-suite-python-ng mcxtrace-suite-python-ng`)

Thank you

Questions?

mcstas.org

mcxtrace.org



SUPPORTED BY



Peter Willendrup · pkwi@fysik.dtu.dk · DTU Physics